


Medium Right Shoulder Sleeve: #510422
Large Right Shoulder Sleeve: #510432
***Right Shoulder Heat Exchanger: #520422**
*fits both Sleeve sizes

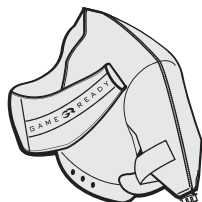

Medium Left Shoulder Sleeve: #510424
Large Left Shoulder Sleeve: #510434
***Left Shoulder Heat Exchanger: #520424**
*fits both Sleeve sizes

Thank you for choosing the Game Ready System. To assure you have the best possible experience please be sure to read the entire User's Manual prior to the use of this product.

Warnings

- Always consult with your physician before using the Game Ready System for any injury. Your physician can determine whether or not the System is appropriate for your injury, and can provide advice on how to best use the System to treat your injury.
- If used incorrectly, cold therapy can cause frostbite or other tissue damage. Always follow the protocol provided by your physician.
- Do not place directly against open wounds, sores, rashes, infections, or stitches. May be applied over clothing or dressing.
- If use of the Game Ready System causes increased pain, swelling, discoloration or loss of sensation, discontinue use and consult with your physician.
- If pressure is uncomfortable, pause unit, and restart with lower pressure.

Care: Gently remove Heat Exchanger from the Sleeve. Hand or machine wash the Sleeve in cold water and mild detergent, or antibacterial soap. Hang to dry. Hand wash the Heat Exchanger, do not machine wash.

Warranty Information:

Sleeve: In case of manufacturer defect, Sleeve may be returned within 7 days of purchase.


Heat Exchanger: 1 year from date of purchase. See warranty card included with Heat Exchanger.




CoolSystems, Inc.
1201 Marina Village Pkwy Suite 200
Alameda, CA 94501
1.888.GAMEREDY

EC REP
EMERGO EUROPE
Molenstraat 15
2513 BH, The Hague
The Netherlands


Shoulder Wrap Application (Right Shoulder Wrap illustrated)

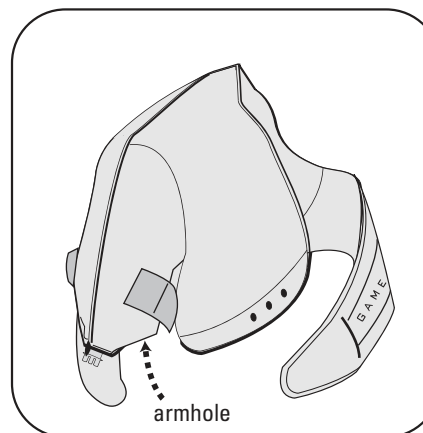


Step 1: Attach the elastic strap to make the armhole.

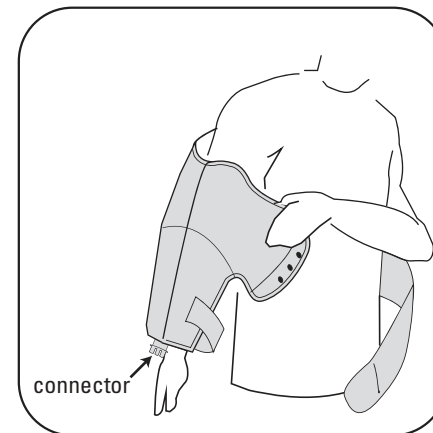

Step 2: Next, slide your arm through the armhole. The connector should hang alongside your arm. You may prefer to connect the hose before putting on the Wrap.

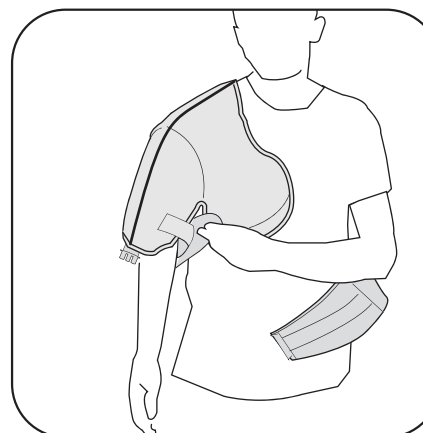

Step 3: Reach under your arm that is in the Wrap, grab the strap and pull it to the front of the Wrap and fasten near underarm on chest.


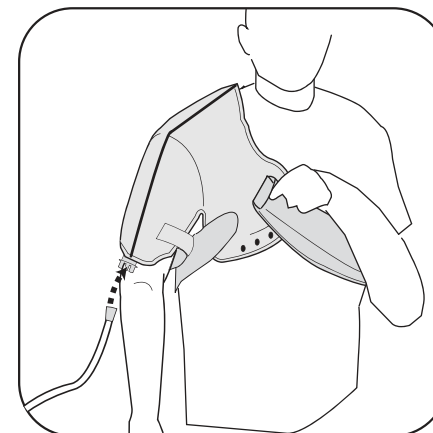


Step 4: Grab the large strap hanging down your back and pull it across your chest on the opposite side of the Wrap and fasten securely to the front of the Wrap. (Larger chest may need the built-in chest extension)

Note: For maximum cooling, place Wrap against bare skin.


Caution: Do not open zipper while Wrap is in use, or applied to the body.

Arm Sling Attachment (included but optional)

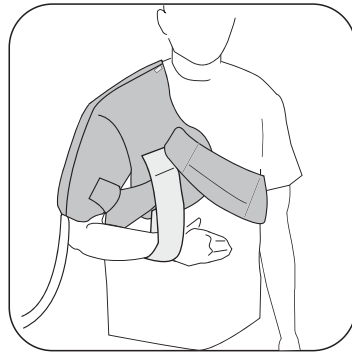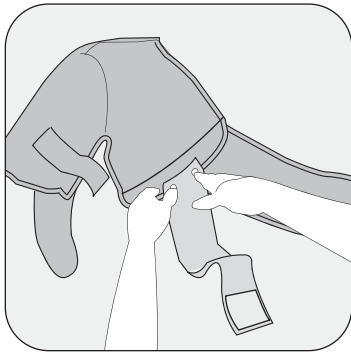

For shoulder injuries that require limited mobility during treatment, an optional elastic sling is included for your use. Simply attach the snaps of the sling at the bottom edge of the front of the Sleeve, and attach the fastener end at a comfortable location on the chest area.

Chest Extension (included but optional)

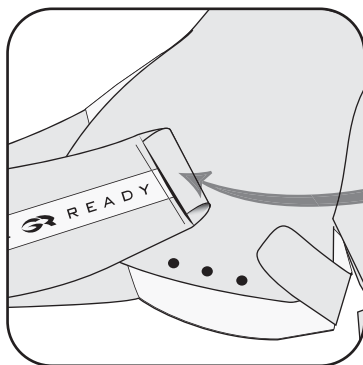

For your convenience, the front chest strap of your Shoulder Sleeve has a built-in chest extension for larger chest sizes. This fold-away extension tucks into a pocket when not needed. The additional panel offers an extra 5 inches of length.

Shoulder Wrap = (Right Shoulder Wrap illustrated)

Shoulder Sleeve

+

Shoulder Heat Exchanger

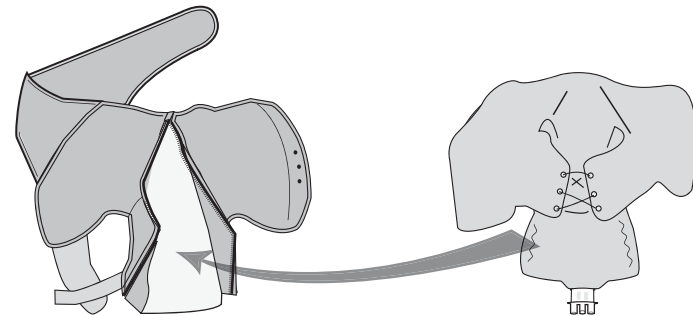


fig. 1 Place Sleeve on flat surface with logo side up and open the zipper. You may prefer to drape it over your knee to insert, as the Wrap has a curve to it.

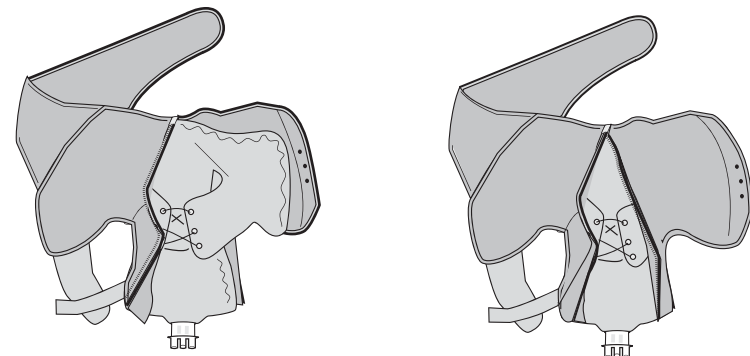


fig. 2 Insert Heat Exchanger into the Sleeve with blue side of the Heat Exchanger down (Blue side of Heat Exchanger should lay against blue side of Sleeve).

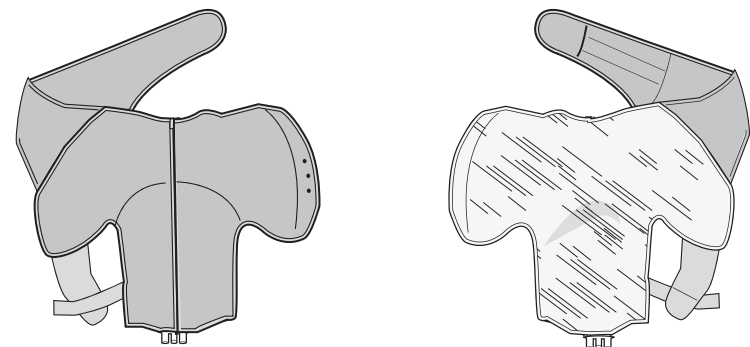


fig. 3 Make sure the Heat Exchanger is flat inside the Sleeve with no folds or creases. Zip up the Sleeve. After assembly, lay the blue side of the Sleeve up, make sure the Wrap is smooth against the inside of the Sleeve.